


# Custom Spinal TLSO Corset Front Orthosis

With Posterior Section and Axilla Straps



Anterior View

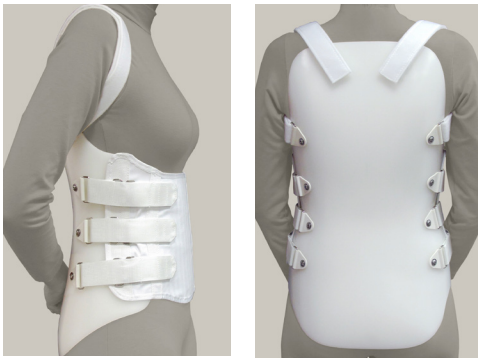

Lateral View

Posterior View

## Product Overview

This TLSO is unique in that it utilizes a soft fabric corset which sits under the chest, supporting the front torso. The back consists of a high profile plastic panel. This combination brace is suitable for multiple diagnosis and is conveniently adjusted. Posterior Section with Axilla Straps (Pictured).

## Spinal Indications

- Herniated Disc
- Low Back Pain
- Compression Fracture
- Lumbar Instability
- Spondylolisthesis
- Osteoporosis
- Degenerative Disc Disease

## Features

- Suitable for patients unable to tolerate a solid plastic anterior panel.
- Adjustable hook and loop closure straps.
- Shoulder straps maintain extension.
- Ease of donning and doffing.

Send your orders to:  
[orders@spinaltech.com](mailto:orders@spinaltech.com)

Send your CAD files to:  
[cad@spinaltech.com](mailto:cad@spinaltech.com)

For product options see  
reverse side.

## **Product Options**

### **Recommended Plastic**

Polyethylene (PE) 3/32", 1/8", 5/32", 3/16" (Least Rigid)

Modified Polyethylene (MPE) 1/8", 5/32", 3/16" (Semi Rigid) (if Transfer pattern is requested)

---

### **Available Plastics**

Copolymer (CP) 1/8" or 5/32" depending on patient size  
(Call for available color options)

---

### **Recommended Aliplast Liner**

3/16" Aliplast

---

### **Other Liner Thickness Available**

1/8", 1/4" Aliplast

Note: This orthosis can also be unlined.

---

### **Spinal Reliefs**

Surgical relief (wide and gradual)

Spinous Process (narrow)

---

### **Fasteners**

Screws and Barrel Nuts

Speedy Rivets

---

### **Other available options**

Axilla straps

Paraspinal reinforcements